



PROCEEDINGS OF THE

27th International Symposium
on Analytical and Environmental Problems

Szeged, Hungary
November 22-23, 2021



University of Szeged

Edited by:
Tünde Alapi
Róbert Berkecz
István Ilisz

Publisher:
University of Szeged, H-6720 Szeged, Dugonics tér 13,
Hungary

ISBN 978-963-306-835-9

2021.
Szeged, Hungary

APPLICATION OF T4 COVALENT FUSION FOR STRUCTURAL STUDIES OF BACTERIAL MEMBRANE TRANSPORTERS

Gerda Szakonyi, Emese Zsuzsanna Galgóczi, Eszter Rajki

University of Szeged, Faculty of Pharmacy, Institute of Pharmaceutical Analysis, H-6720

Szeged, Somogyi u. 4.

e-mail: szakonyi.gerda@szte.hu

Investigation of the transmembrane proteins found in bacteria is important, among others, to gain knowledge how to effectively control certain pathogens. Molecular chaperones such as T4 lysozyme or antibodies are found to be efficient for studying the structure of membrane proteins. Our aim is to covalently fuse T4 lysozyme to the C terminus of *Hemophilus influenzae* HI0610 membrane protein to study the structure of this fusion protein by X-ray crystallography after cloning, expression, purification and crystallization.

HI0610 transmembrane protein was cloned into pTTQ18 plasmid by using recombinant DNA technology. T4 lysozyme was covalently fused to the target protein during the cloning. The protein expression was followed by SDS gel-electrophoresis to investigate the presence of the required macromolecule. Optimization of protein expression was also performed where IPTG concentration, temperature and time dependence were tested. Subsequently, the protein was purified by affinity chromatography and the efficiency of Ni-NTA and Talon resins was compared. The target protein was further purified by size exclusion and ion exchange chromatography.


T4 lysozyme was successfully fused to HI0610 fucose-proton symporter by genetic engineering. Gel-electrophoresis confirmed that the target macromolecule was expressed in an appropriate amount. Optimal conditions for protein expression were the following in Luria-Broth medium: 37°C, 12 h of induction time, 0.5 mM IPTG. During affinity chromatography, the target protein was recovered in high purity. Investigating the eluted proteins obtained during protein purification using two affinity columns, we concluded that the Ni-NTA resin bound the target protein more efficiently than that of Co²⁺ ion containing Talon resin. However, after affinity chromatography further purification steps were necessary to reach the highest protein purity. Therefore, size exclusion and cation exchange chromatography were also performed to remove the contaminating proteins.

The final analysis of the purified concentrated sample revealed no detectable contaminations even by the sensitive silver staining and become suitable for crystallization and structural studies.

APPLICATION OF T4 LYSOZYME COVALENT FUSION FOR STRUCTURAL STUDIES OF MEMBRANE TRANSPORTERS

Gerda Szakonyi, Emese Zsuzsanna Galgóczi, Eszter Rajki

University of Szeged, Faculty of Pharmacy
Institute of Pharmaceutical Analysis

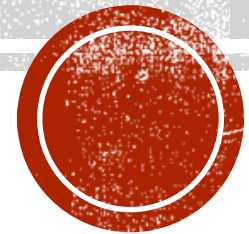



INTRODUCTION

- Aim: structural and functional characterization of *Hemophilus influenzae* membrane transporter
- L-fucose/proton symporter: **HI0610**
- Major Facilitator Superfamily (MFS)
 - Homology with eukaryotic proteins
 - Potential pharmaceutical drug target

Successful expression



Failed crystallization

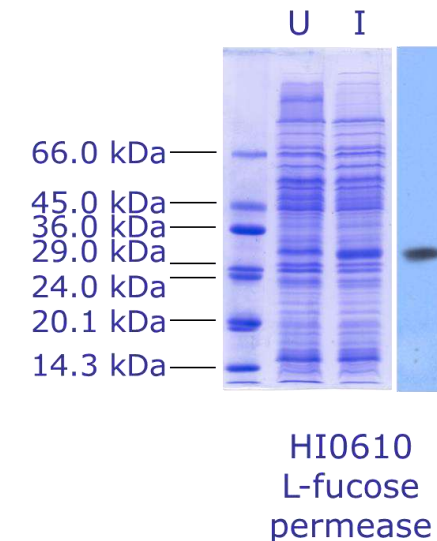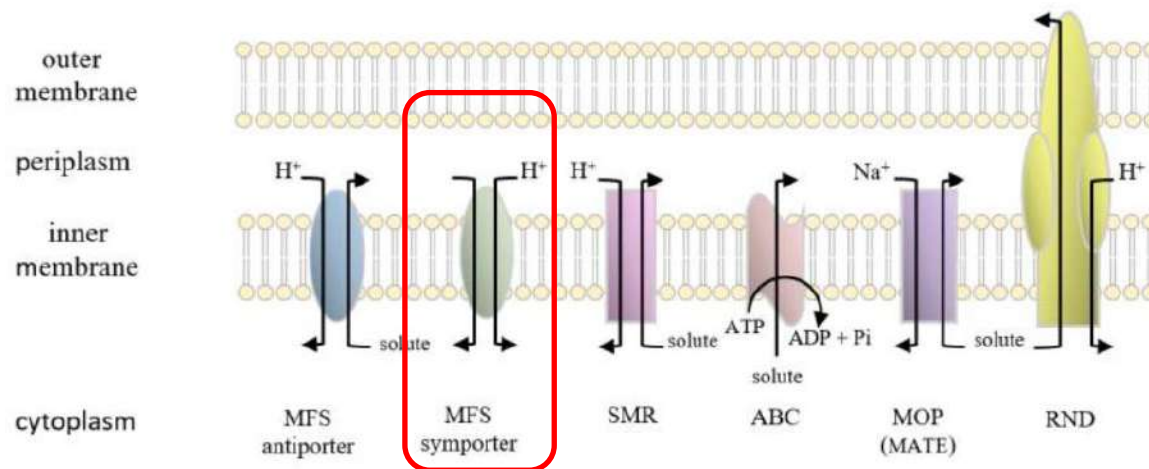


STRUCTURAL DETERMINATION OF MEMBRANE PROTEINS

- Difficulty in membrane protein crystallization:
 - Limited possibilities for protein protein interactions → increasing the structurally ordered extramembrane domains
- Solution: crystallization chaperon protein
 - covalent fusion → T4 lysozyme

Protein	Chaperone	PDB code	Resolution (Å)
A_{2A} Adenosine receptor	T4 lysozyme	3EML	2.6
A₂ AR	T4 lysozyme	3QAK	2.7
β₂ Adrenergic receptor	T4 lysozyme	2RH1	2.4
β₂AR-G_s complex	T4 lysozyme	3SN6	3.2
CXCR chemokine receptor	T4 lysozyme	3ODU	2.5
Dopamine D3 receptor	T4 lysozyme	3PBL	2.9
Histamine H₁ receptor	T4 lysozyme	3RZE	3.1
β₂AR-G_s complex	T4 lysozyme (non-covalent)	3SN6	3.2

- Topology prediction:
 - SOSUI program (Nagoya University, Japan)
 - Affinity group: His-Tag

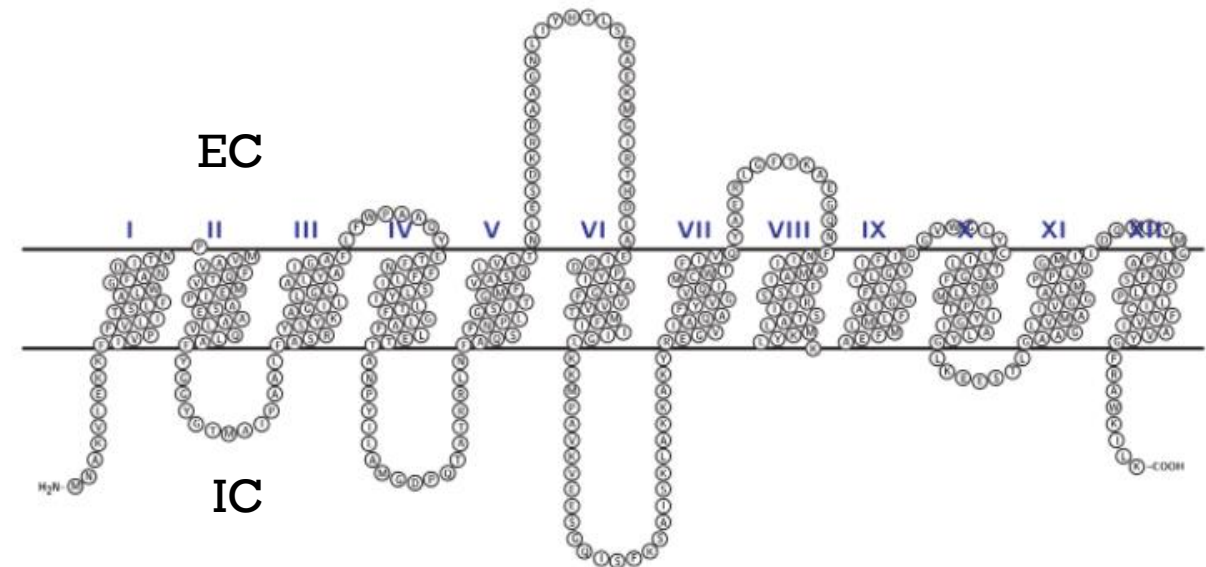

METHODS

- Recombinant pTTQ18 – HI0610 – T4



Transformation (heat shock)

- BL21(DE3) *E. coli* strain



IPTG: Isopropyl β -D-1-thiogalactopyranoside

- Expression




Water lysis


- Purification

Long term goal:

- High scale purification
- Functional characterization
- Crystallization
- Structural determination



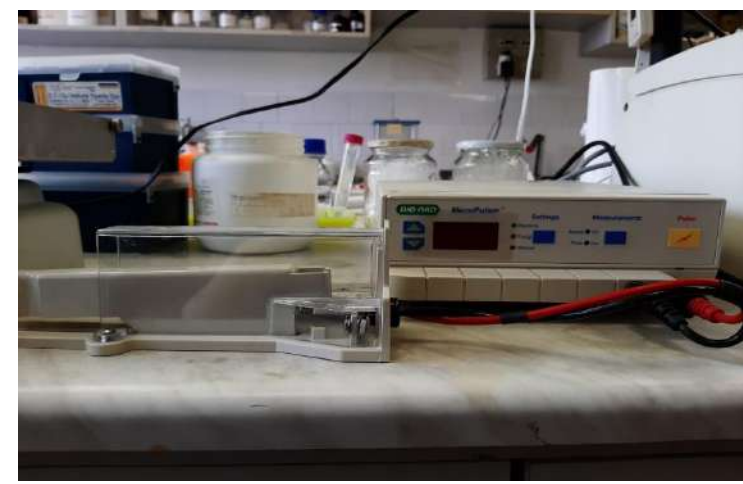


Electroporation

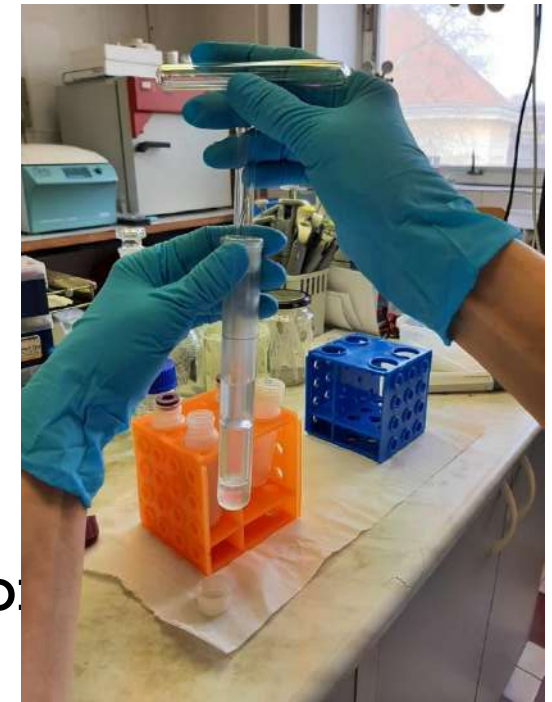


Water lysis

GENETIC ENGINEERING

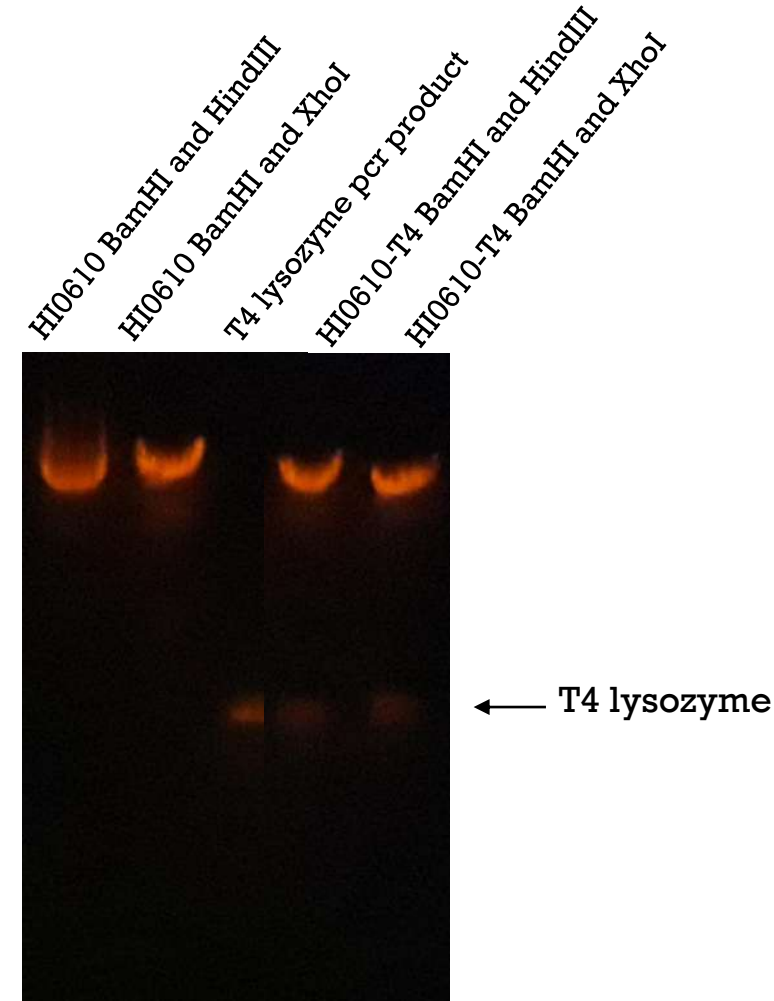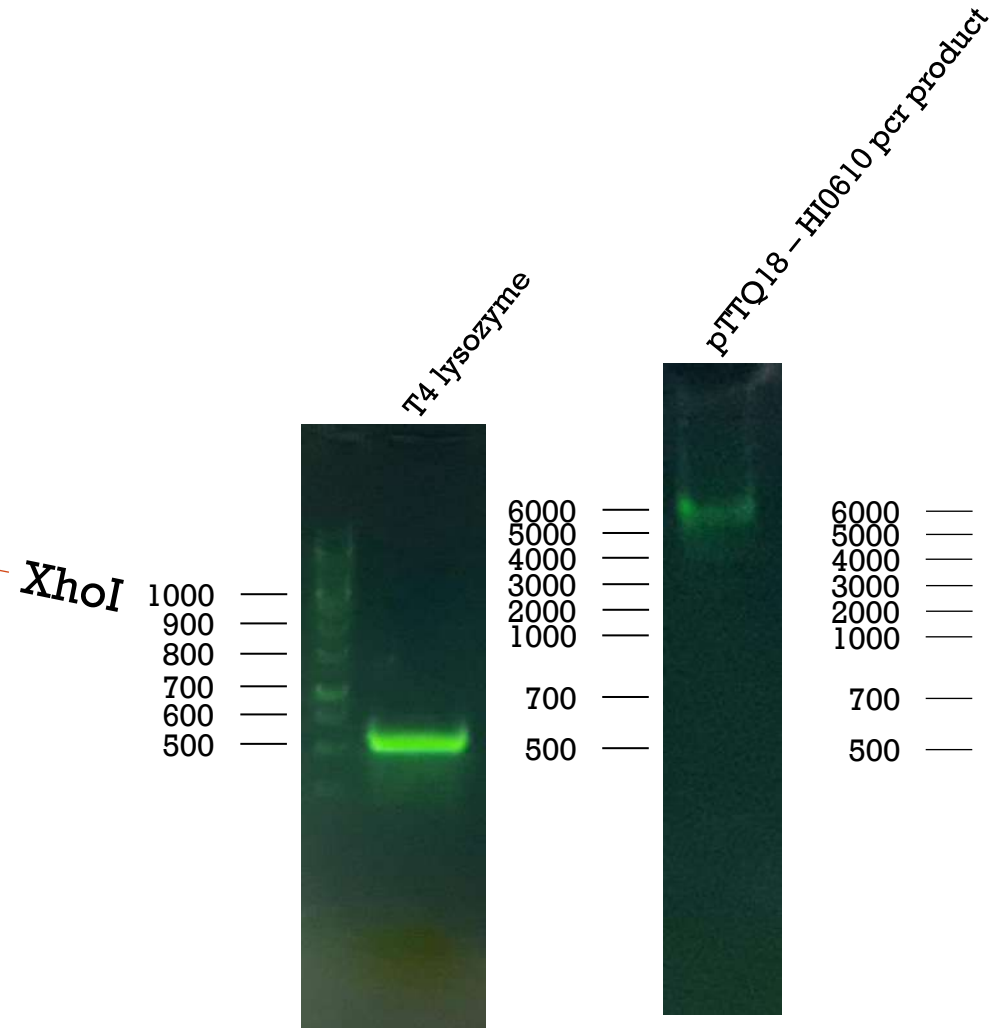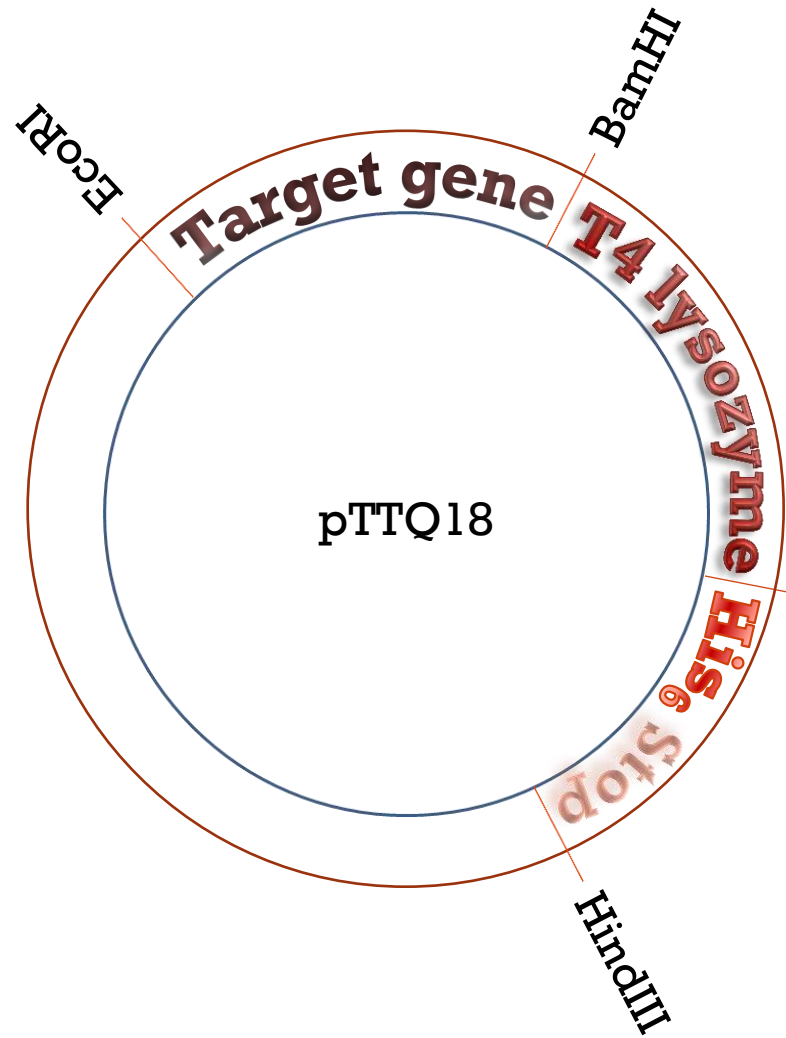


RESULTS – SDS-PAGE

○ Expression test

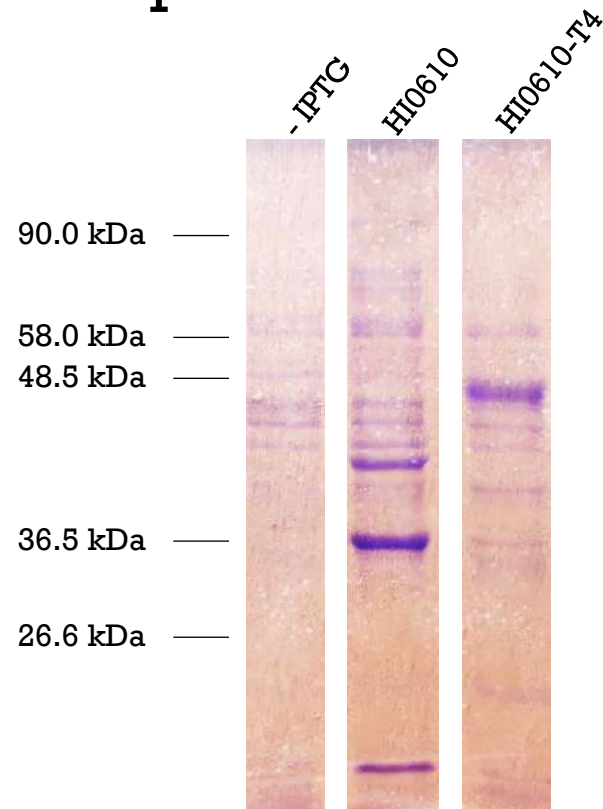

- LB media
- 37 °C
- 3 h induction


○ Optimization of IPTG concentration

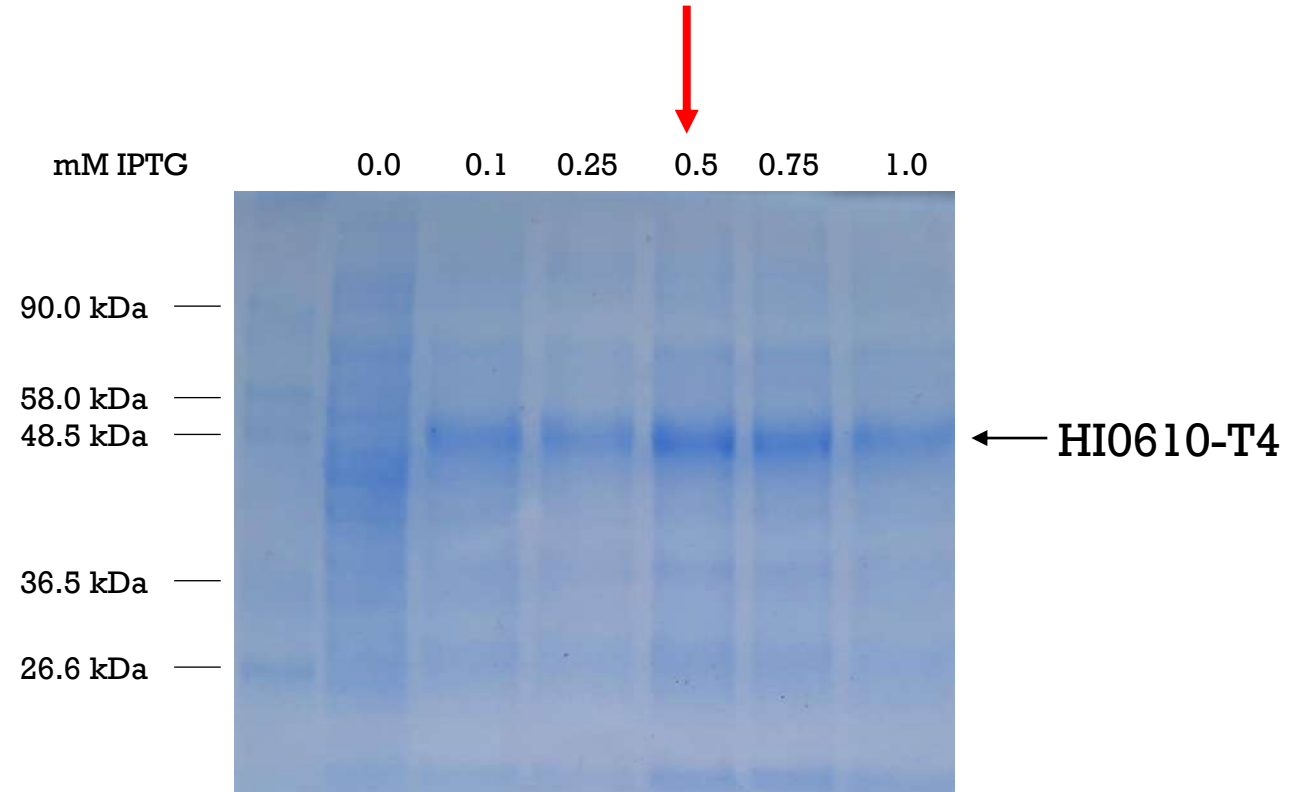

- Optimal IPTG concentration: 0.5 mM


OPTIMIZATION OF INDUCTION TIME

○ $T=25^{\circ}\text{C}$

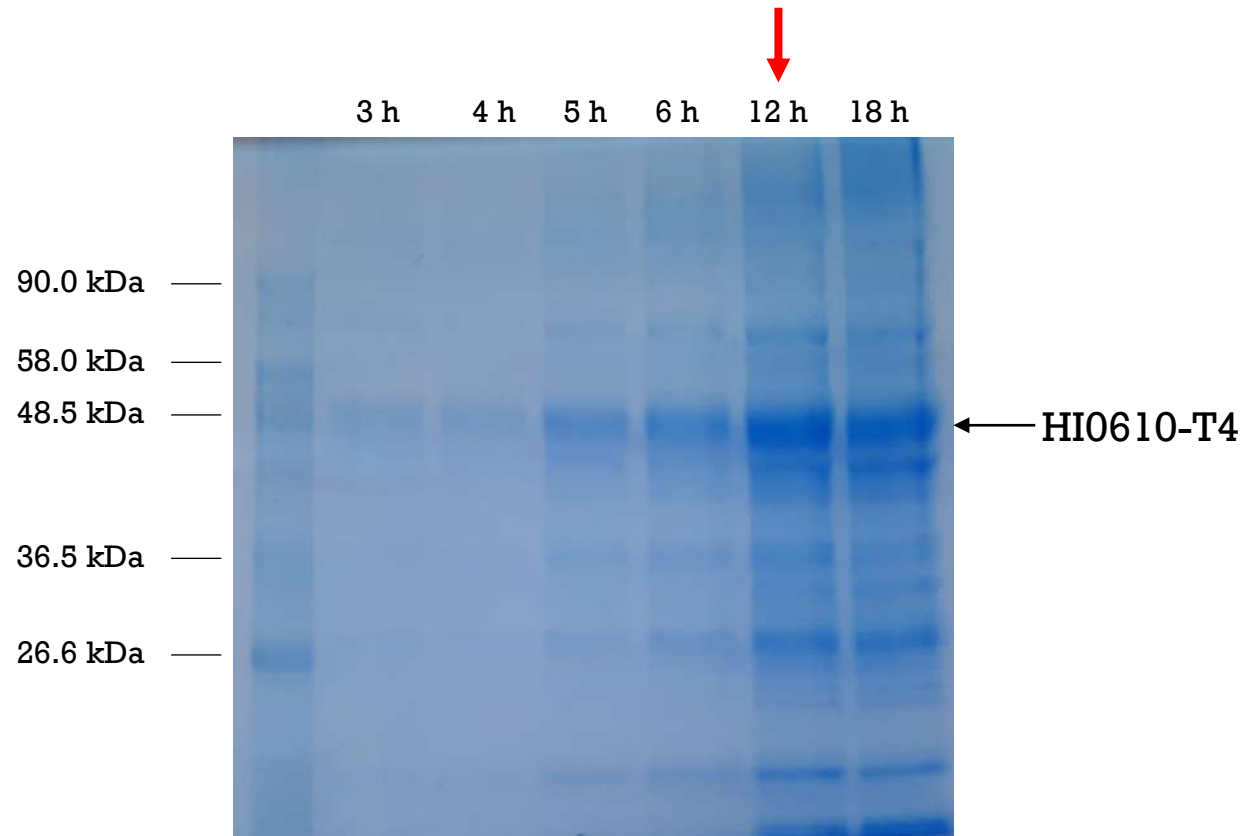

○ Optimal induction time: 12 h


○ $T=37^{\circ}\text{C}$

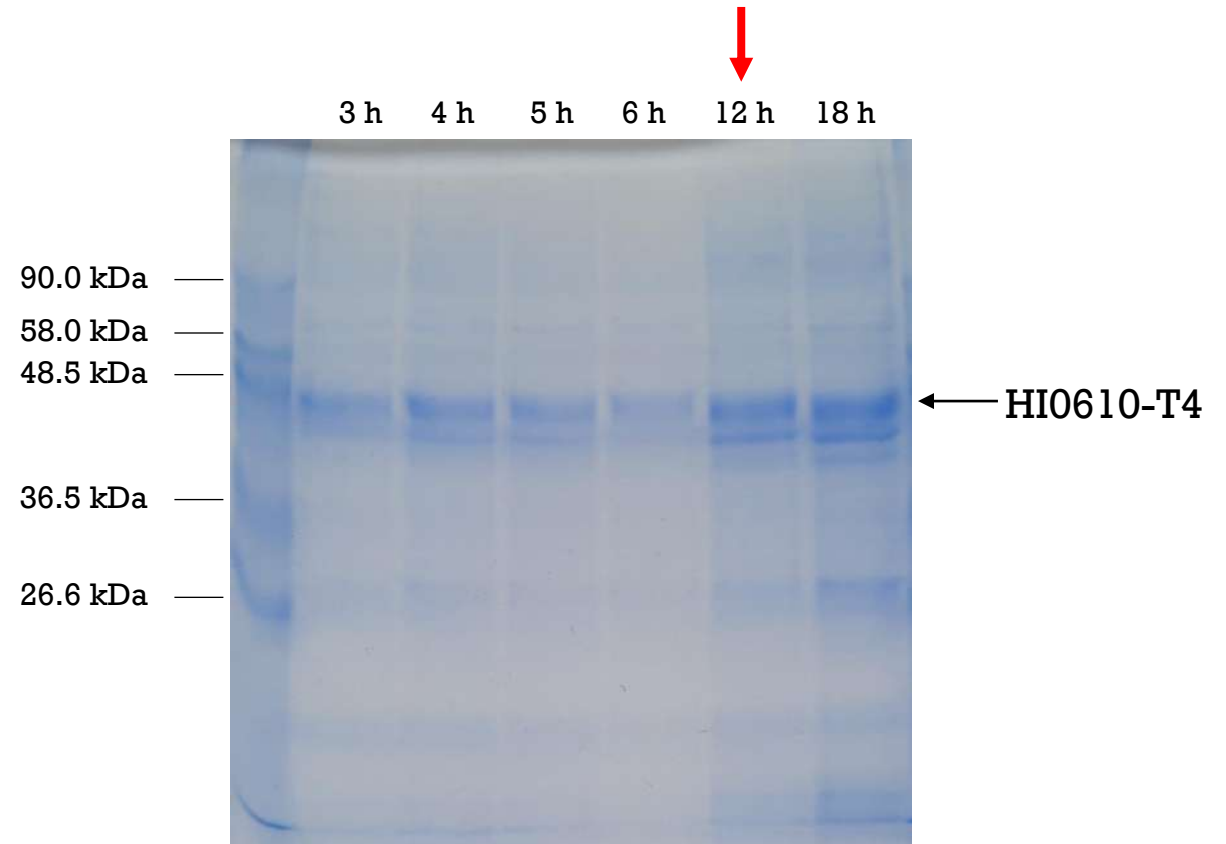


○ Optimal induction time: 12 h

AFFINITY CHROMATOGRAPHY

- Purification: affinity chromatography

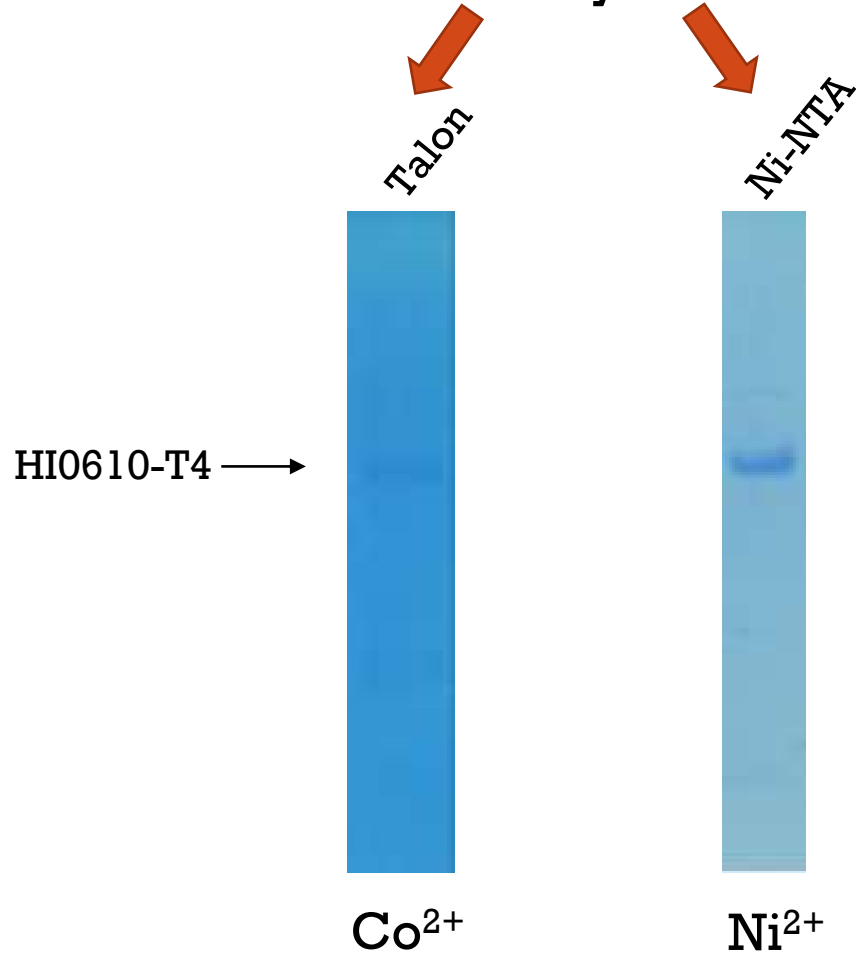

- Sample after Ni-NTA affinity chromatography


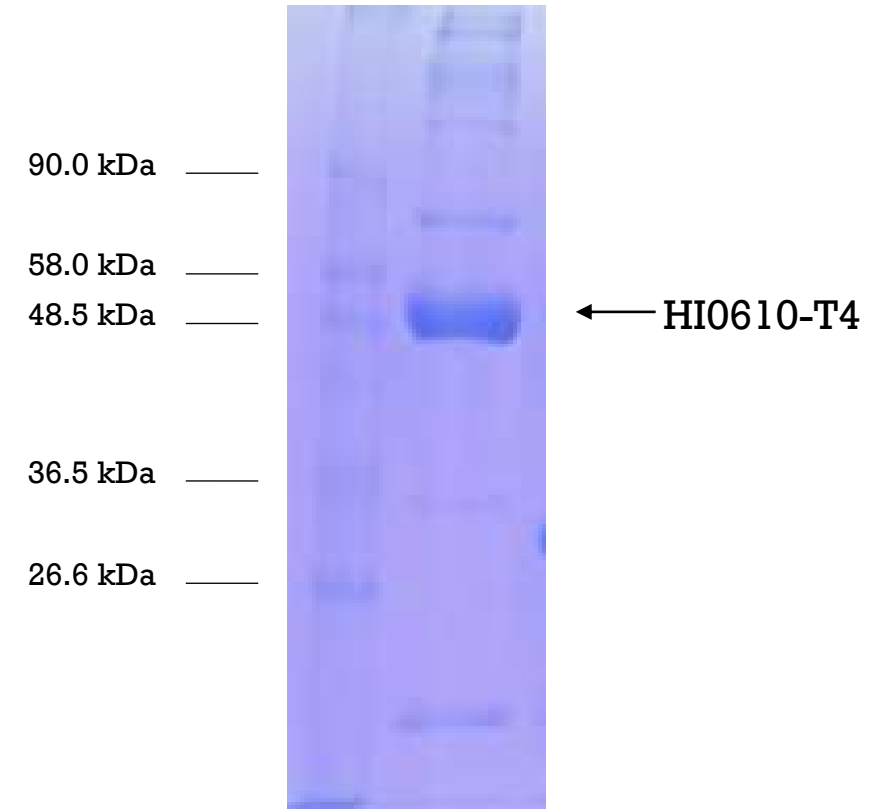


SIZE EXCLUSION CHROMATOGRAPHY

- ÄKTA FPLC system, Superdex 200 column

- Fractions on SDS PAGE



20 mM Tris pH 8,0
150 mM NaCl
5% glycerol
0.05% DDM

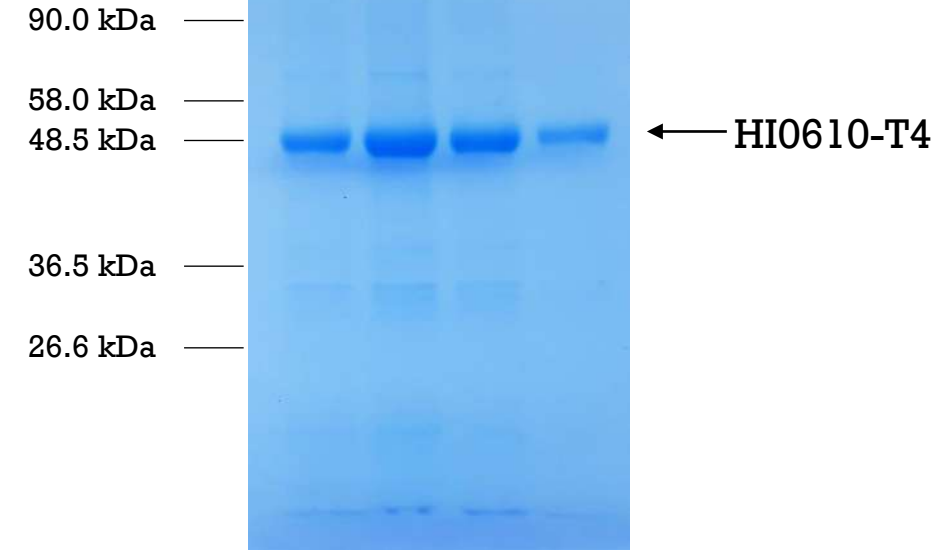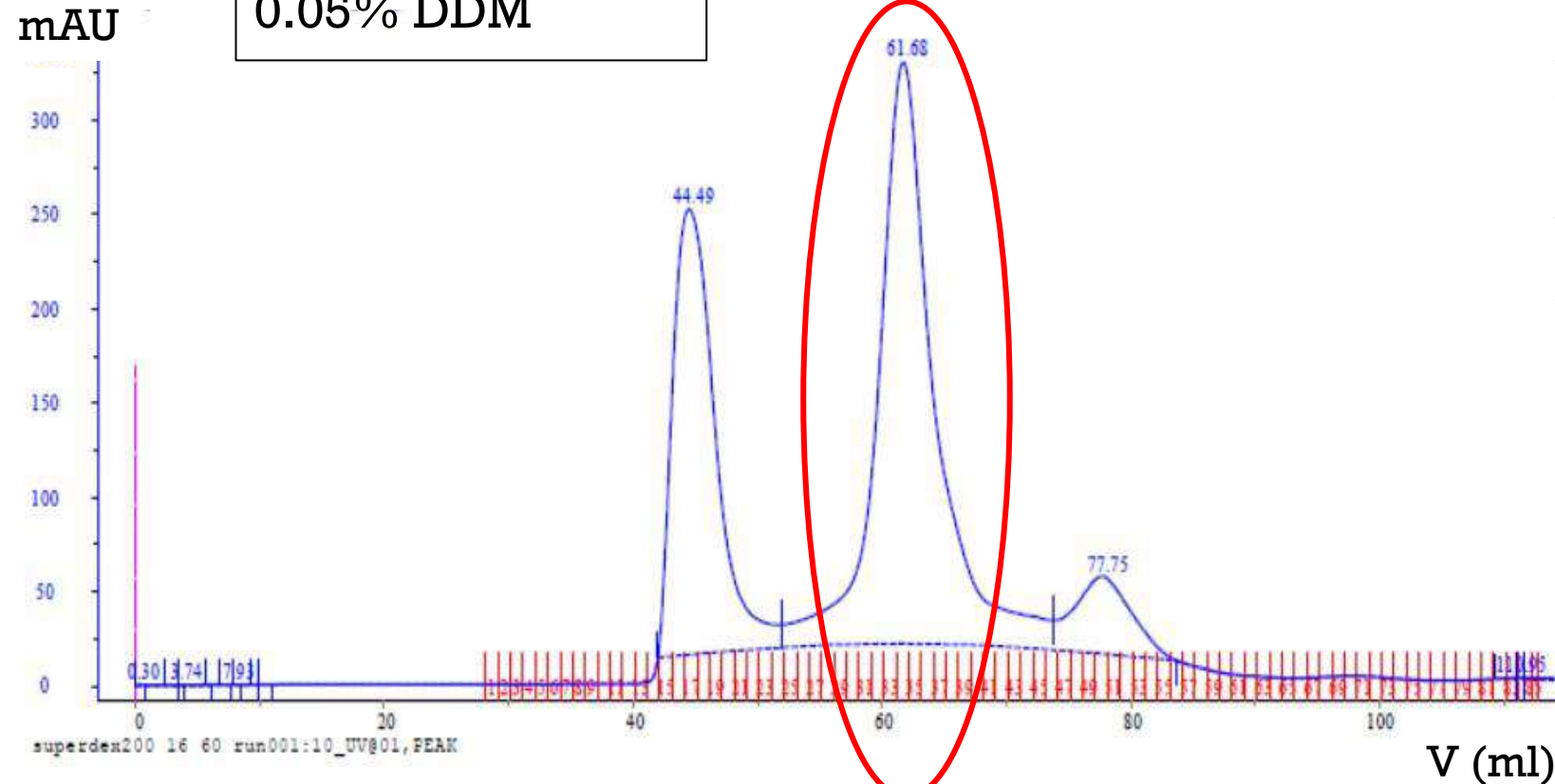

ION EXCHANGE CHROMATOGRAPHY


- ÄKTA FPLC system, Resource S column

- Target protein – SDS-PAGE

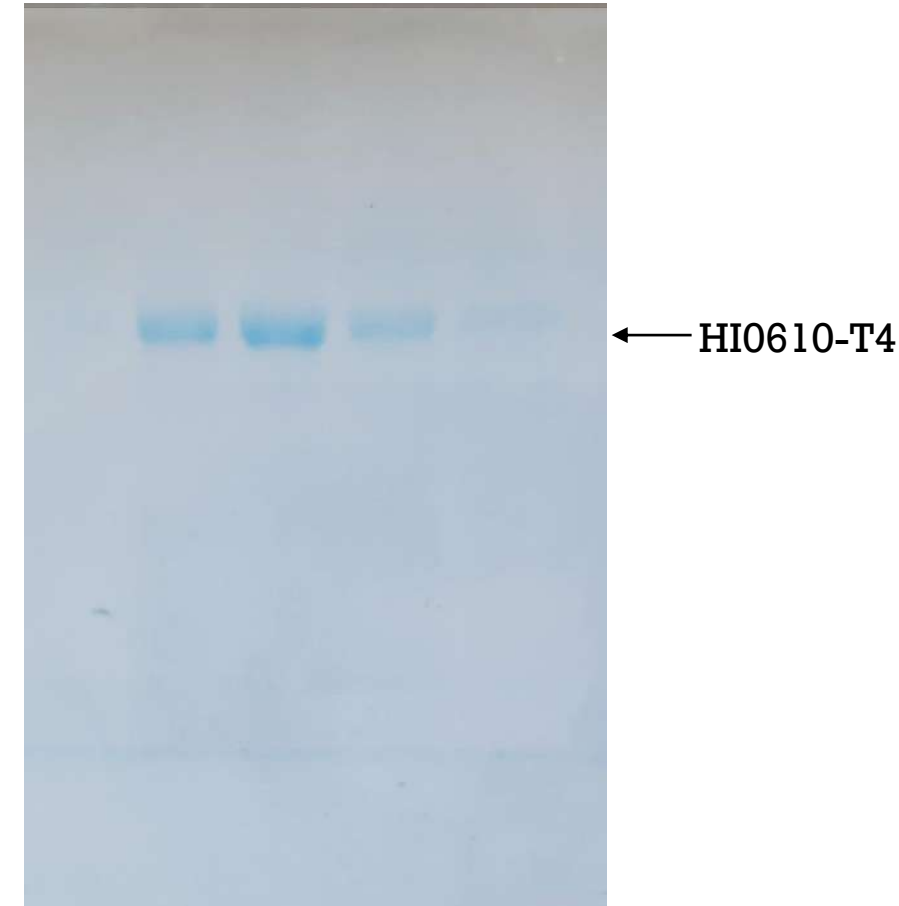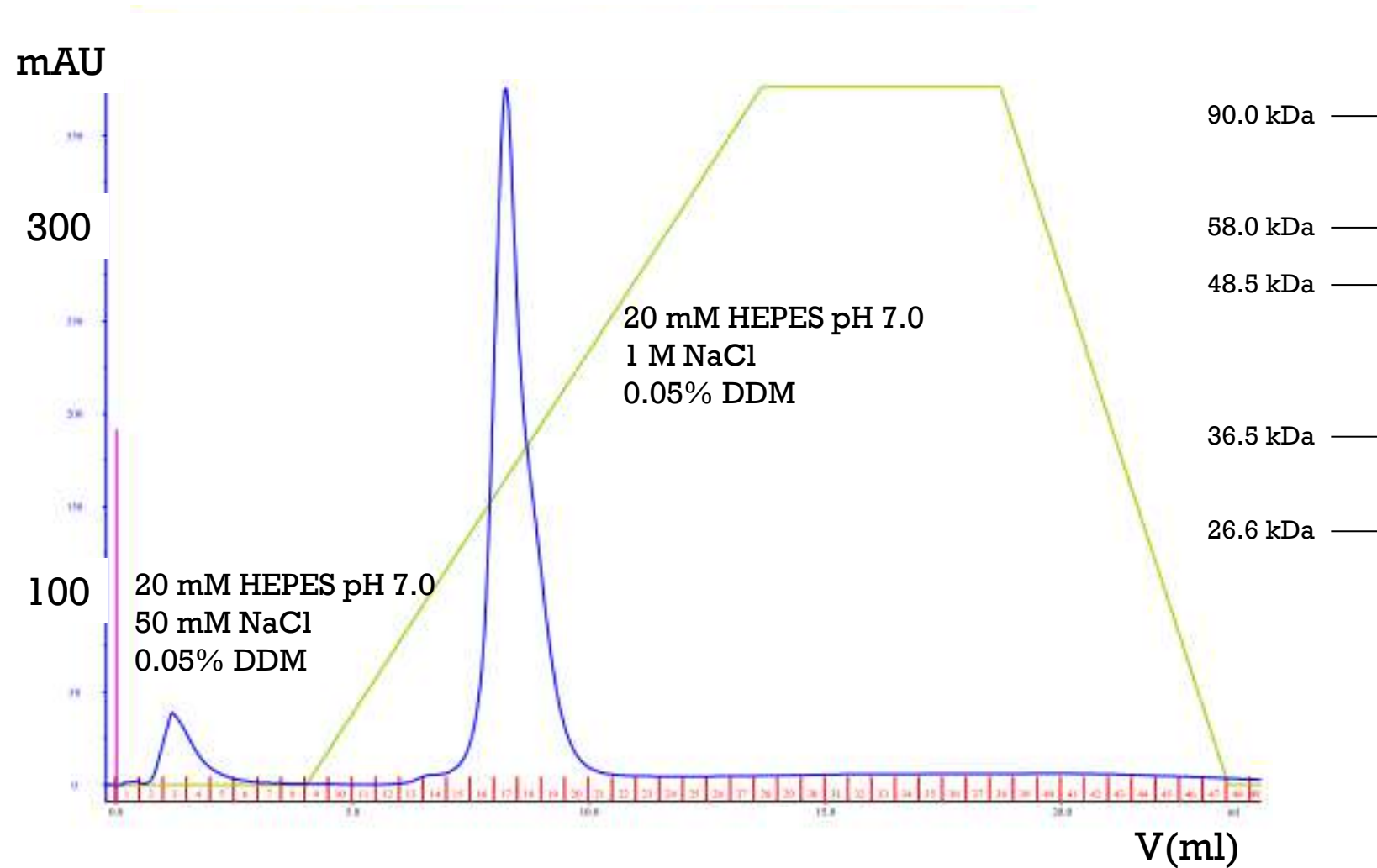


CONCENTRATED SAMPLE ON SDS-PAGE

Coomassie Brilliant Blue

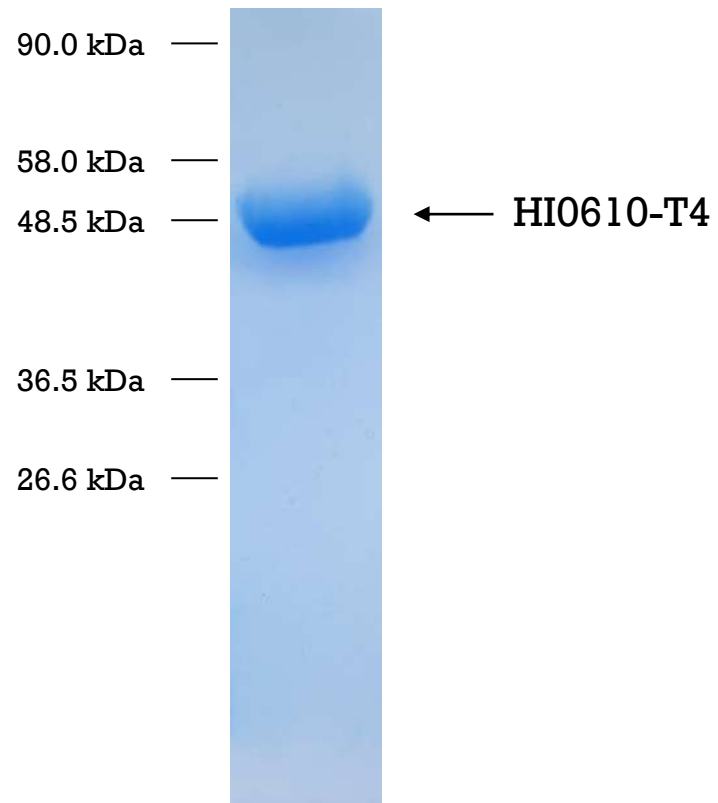

Silver staining

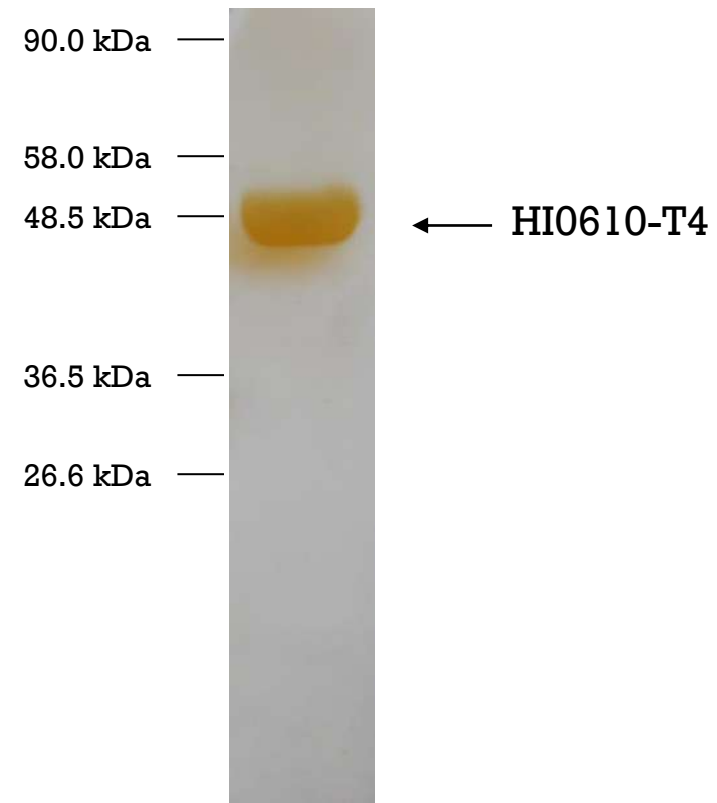

DISCUSSION

- Successful HI 0610 target protein C-terminal T4 lysozyme fusion
- The optimal conditions of protein expression and purification are:
 - 0.5 mM IPTG
 - 37°C
 - LB medium
 - 12 h
 - Ni²⁺ resin
- The application of size exclusion and ion exchange chromatography resulted in high purity of protein

ACKNOWLEDGEMENT

Lóránt Lakatos D. Sc. (BRC)

Erika Kunosné Tóth

István Ilisz D. Sc.